



# **Laboratory Manual for Aboveground Understory and Soil Organic Carbon Analysis**

Department of Forest & Park Services  
Ministry of Agriculture & Forests





# Table of Contents

|   |          |
|---|----------|
| FOREWORD .....  | 2        |
| ACRONYMS AND ABBREVIATIONS .....  | 4        |
| <b>Chapter I:</b>   |          |
| <b>Introduction .....</b>   | <b>1</b> |
| <b>Chapter II:</b>  |          |
| <b>Procedures for aboveground understory plant carbon assesment .....</b> | <b>2</b> |
| 1. Samples sorting by client.....   | 2        |
| 2. Samples reception and registration by the Laboratory staff.....        | 3        |
| 3. Sample processing: Oven drying and dry matter estimation.....          | 4        |
| 4. Estimating carbon content in aboveground understory.....               | 5        |
| <b>Chapter III:</b>   |          |
| <b>Procedures for soil organic carbon analysis .....</b>                  | <b>7</b> |
| 1. Soil sample sorting by the client .....                                | 7        |
| 2. Samples reception and registration by Laboratory staff.....            | 7        |
| 3. Soil samples processing.....   | 8        |
| 4. Soil moisture content determination and correction factor.....         | 10       |
| 5. Analysis for soil Organic Carbon using Walkley and Black Method.....   | 12       |
| 6. Analysis of Soil Organic Carbon using CHN628 Analyzer.....             | 16       |
| <b>Figures</b>  |          |
| Fig.1.Samples in a sack .....   | 2        |
| Fig.2a Samples sorting.....   | 2        |
| Fig.2b. Samples arranged.....   | 3        |
| Fig.3. Filled sample submission Form .....                                | 3        |
| Fig.4. Plant Sample Register.....   | 3        |
| Fig.5. SPAL's Database computer .....                                     | 4        |
| Fig.6. Electronic balance.....  | 4        |
| Fig.7. Weight of an empty tray.....                                       | 4        |
| Fig.8. Tray plus fresh sample.....  | 4        |
| Fig.9. Tray plus fresh sample in oven.....                                | 4        |
| Fig.10.Weight of tray plus dry sample .....                               | 5        |
| Fig.11 a. Record weight of dry sample.....                                | 5        |
| Fig.11b. Complete weight record.....                                      | 5        |
| Fig.12. Bulk soil samples.....  | 7        |

|  |    |
|--|----|
| Fig.13. Sorting sample.....                            | 7  |
| Fig.14 Sample Submission Form .....                    | 7  |
| Fig.15. Soil Registration Log book .....               | 8  |
| Fig.16. SPAL's Database computer.....                  | 8  |
| Fig.17. Samples for drying.....                        | 8  |
| Fig.18. Drying soil samples .....                      | 8  |
| Fig.19. Stones Seggregated from soil and weighed ..... | 9  |
| Fig.20a. Sieving soil sample.....                      | 9  |
| Fig.20b. Grind sample with mortar and pestle.....      | 9  |
| Fig.21. Taking soil weight .....                       | 9  |
| Fig.22. Stone & soil weight recorded.....              | 9  |
| Fig.23. Soil sample bottled .....                      | 9  |
| Fig.24. Soil sample bottled with code .....            | 9  |
| Fig.25. Printed soil Carbon results.....               | 9  |
| Fig.26 Weighing for moisture .....                     | 11 |
| Fig.27. Oven drying for moisture.....                  | 11 |
| Fig.28. Oven dried sample.....                         | 11 |
| Fig.29. Cooling in desiccators .....                   | 11 |
| Fig.30. Weighing dry sample.....                       | 11 |
| Fig.31. Dry weight recorded.....                       | 11 |
| Fig.32. Moisture calculation in Database .....         | 11 |
| Fig.33. Weighing for carbon analysis.....              | 14 |
| Fig.34. Adding sulphuric acid to sample.....           | 14 |
| Fig.35. Titration .....                                | 14 |
| Fig.36. Record the reading .....                       | 14 |
| Fig.37. Result calculation in Database .....           | 14 |
| Fig.38. CHN628 Analyzer .....                          | 16 |
| Fig.39. PC to control the system .....                 | 16 |
| Fig.40 Gases cylinders .....                           | 16 |
| Fig.41. Balance with sample cup.....                   | 17 |
| Fig.42. The carousel.....                              | 17 |
| Fig.42a. Loading carousel .....                        | 17 |
| Fig.43. Printed result .....                           | 17 |

## Tables

|              |    |
|--------------|----|
| Table 1..... | 5  |
| Table 2..... | 5  |
| Table 3..... | 6  |
| Table 4..... | 15 |
| Table 5..... | 18 |

# FOREWORD

I am very happy that the collaboration between our department and the Department of Agriculture has resulted in the development of this laboratory manual for organic carbon analysis of plants and soil.

This laboratory manual fulfills three important purposes:

1. It serves as an important reference for organic carbon analysis of both plant and soil in the country
2. It provides institutional support for pursuing periodic national forest inventory
3. It expands the scientific and technical capacity of both forestry and agriculture sector

This is one of the important building blocks for institutionalizing **measuring, reporting and verification** (MRV) of forest carbon and consequently fulfilling our commitment to remain carbon neutral.

For this, I would like to thank the Department of Agriculture, and particularly Dr. Karma Dema, Program Director of National Soil Service Centre and Mr. Jamyang, Specialist at Soil and Plant Analytical Laboratory for being receptive to the collaboration. I remain highly appreciative and grateful for the support and hard work rendered by Mr. Jamyang and his team at SPAL in analyzing almost 2000 samples of plant and soil collected from our NFI plots and for bringing out this laboratory manual. The forest carbon assessment wouldn't be complete without the inclusion of above ground understory and soil carbon and therefore, their contribution to the effort is very significant.

I would like to congratulate Mr. Lobzang Dorji, Chief Forestry Officer at Forest Resources Management Division (FRMD) and his NFI team for successfully leading the collaboration with NSSC and for jointly producing this manual. I understand this manual highlights the expansion of forestry sector beyond timber and timber volume;and enhances our capacity to account for the role of forests in sequestering carbon.

I would also like to share my appreciation to BMUB, GIZ and ICIMOD for supporting our efforts on forest carbon assessment through the REDD+ Himalayas project and for supporting the development of such knowledge products.

This knowledge product will be useful to many other agencies, students, scholars and researchers interested in above ground understory plant and soil carbon analysis procedures.

A handwritten signature in blue ink, appearing to read 'Phento Tshering', with a long horizontal stroke extending to the right.

Phento Tshering  
Director

## ACRONYMS AND ABBREVIATIONS

|                |  |
|----------------|--|
| ABG            | Aboveground  |
| BMUB           | Federal Ministry for the Environment, Nature Conservation, Building and Nuclear Safety                     |
| BT FEC         | Bhutan Trust Fund for Environmental Conservation   |
| CI             | Confidence Interval  |
| cm             | Centimeter   |
| °C             | Degree centigrade/celcius  |
| DoFPS          | Department of Forests and Park Services  |
| EU             | European Union   |
| FAO            | Food and Agriculture Organization  |
| FCPF           | Forest Carbon Partnership Facility   |
| FRMD           | Forest Resources Management Division   |
| g              | Gram   |
| GCCA           | Global Climate Change Alliance   |
| GIS            | Geographic Information System  |
| GIZ            | Deutsche Gesellschaft für Internationale Zusammenarbeit - German Corporation for International Cooperation |
| GS             | Growing Stock  |
| Ha             | Hectare (unit of area)   |
| ICIMOD         | International Centre for Integrated Mountain Development   |
| JID            | Job Identification Number  |
| KG or kg       | Kilogram   |
| KM or km       | Kilometer  |
| m <sup>3</sup> | Cubic meter (unit of volume measurement)   |
| Mc             | Moisture correction factor   |
| ml             | Milli liter  |
| Mg/ ha         | Megagrams per hectare  |
| Mg C/ha        | Megagrams of Carbon per hectare or Tonnes of carbon per hectare  |
| MRV            | Measurement, Reporting and Verification  |
| MoAF           | Ministry of Agriculture and Forests  |
| MoE            | Margin of Error  |
| NFI            | National Forest inventory  |
| NFMS           | National Forest Monitoring System  |
| NSSC           | National Soil Service Centre   |
| REDD           | Reducing Emission from Deforestation and Degradation of forests  |
| RGoB           | Royal Government of Bhutan   |
| RNR SSP        | Renewable Natural Resources Sector Support Programme   |
| R-PP           | REDD Readiness Project   |
| SPAL           | Soil and Plant Analytical Laboratory   |
| SNV            | Netherlands Development Agency   |
| UN REDD -TS    | UN REDD Targeted Support   |
| UN             | United Nations   |
| WB             | World Bank   |





# Chapter I: Introduction

**D**eveloped with the intention to provide transparency to the procedures adopted for National Forest Inventory (NFI) of Bhutan (2012-2015) and to serve as institutional reference for future NFI, this is a complementary document to the field manual “A Field Guide for Aboveground Understory (ABG) and Soil Carbon Assessment” that was developed by Forest Resources Management Division in 2014.

This manual provides step by step procedure for handling and processing of plant and soil samples in the laboratory. Together with the field manual, it serves as the complete guidance on carrying out the organic carbon assessment of aboveground understory and soil samples.

This is an important result of the on-going collaboration between Forest Resources Management Division, Department of Forests and Park Services and Soil and Plant Analytical Laboratory (SPAL), National Soil Service Centre, Department of Agriculture of Ministry of Agriculture and Forests (MoAF) on ABG and soil carbon analysis. This manual thus presents the standard analysis procedures being adopted and implemented at SPAL, applied to the overall ABG and soil carbon assessment methodology developed by FRMD and soil scientists of USDA-Forest Service and the University of Washington of USA.

The document will be useful to any future NFI work, academician, researchers and students in the country, interested in studying or conducting a research on organic carbon.

# Chapter II:

## Procedures for aboveground understory plant carbon assesment

### PRINCIPLE:

With climate change and gobal warming issues, it is felt increasingly important to estimate the carbon stock of forest soil and vegetation. This chatppter therefore, deals in assessing the carbon content of aboveground understory plant samples. For this, fresh plant samples need to be processed in the laboratory and dry matter estimated. Plant dry matter refers to the weight of dry sample after removing all moisture at 70°C and is expressed as percentage of fresh weight. An accurate assesment of dry matter is important to measure biomass per specified plot or area, which is essential for estimation of carbon content in aboveground understory sample using a factor of 1.75. In general, approximately half of plant biomass is considered to be made of carbon molecules.

The standard procedure for assesing the dry weight and finally estimating the carbon tonnes per hectare are presented below.



Fig.1.Samples in a sack



Fig.2a. Samples sorting

### 1. Samples sorting by client

- 1.1 Upon reaching plant samples in the laboratory, a client should first do the sorting of the samples which are brought in the sacks Figure 1, and arrange them in sequence (Fig.2a) according to their identification code or labels as listed (Fig.2b) in the Sample Submission Form as

in Figure 3. If not already filled in the Sample Submission Form the client has to fill-in the FORM as he or she wants the results to appear.

If needed SPAL technician can provide helping hand and guide how samples are to be arranged.

- 1.2 After the samples are arranged the client can then hand over them to the Laboratory Sample Receptionist who will check if they are carefully packed in a linen or a plastic bag and clearly labeled for identification. Any sample not properly labeled, or labeled but not clear and discrepancies arise, such samples will be rejected unless the person handing in could clarify the confusion.

The Sample Receptionist will also cross check whether the number of samples and their identification Code listed in the Sample Submission Form tally with the physical count and actual code on the samples

## 2. Samples reception and registration by the Laboratory staff

The Sample Receptionist should cross check whether the number of samples submitted matches with their identification Code listed in the Sample Submission Form (Fig 3) with the physical count and actual code provided in the samples list. Any mismatch should be either instantly corrected or discarded. The Sample receptionist should also check if the contact address provided in the Form is correct

The sample receiving staff should then transcript the information from the Sample Submission Form Fig 3 into the Sample Registration Form Fig 4. This in turn should be entered into the SPAL's Database computer Fig 5 to assign Lab

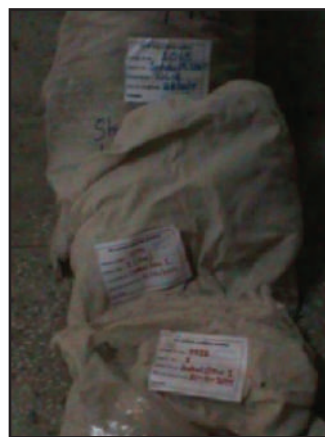


Fig.2b. Samples arranged

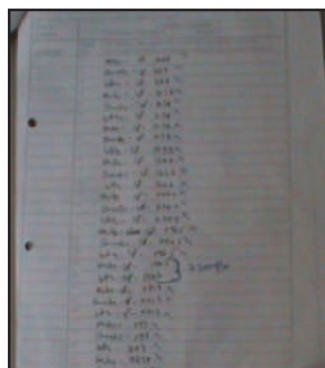


Fig.3. Filled sample submission Form



Fig.4. Plant Sample Register

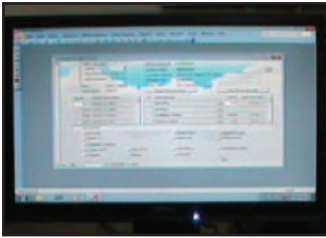


Fig.5. SPAL's Database computer



Fig.6. Electronic balance



Fig.7. Weight of an empty tray



Fig.8. Tray plus fresh sample

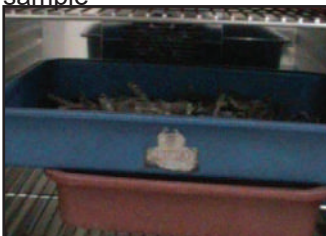


Fig.9. Tray plus fresh sample in oven

Number and Job Identification Number (JID). This JID may be used by the client for future reference for any inquiry. Once lab number is assigned to respective samples the samples are then processed as in section 3 below.

The completed Sample Submission Form is signed and photocopied. One copy is filed and the other is returned to the originator. This is to facilitate the better communication.

### 3. Sample processing: Oven drying and dry matter estimation

Samples received for Carbon analysis are often heterogeneous due to variation in stem, leaf content and the chopping method used. The laboratory attempts to minimize this variation by thorough chopping of larger sample pieces down to a uniform size where possible. In order to estimate the plant Carbon follow the below steps.

- 3.1 Using digital electronic balance (Fig.6), weigh the plastic tray or the linen bag. (Fig.7)
- 3.2 Transfer the sample into the plastic tray or in the linen bag.
- 3.3 Weigh the fresh sample + the container Fig.8.
- 3.4 Record the weight of the wet sample.
- 3.5 Weigh a sub-sample of the wet sample+ container
- 3.6 Put the sample in the ventilated electric oven Fig.9, for drying
- 3.7 Set the temperature to 70 °C and heat over night
- 3.8 Next morning switch off the oven and let it cool.
- 3.9 Weigh the sample + container to obtain the weight of the dry sample Fig 10.
- 3.10 Record all the values as in Fig 11.

#### 4. Estimating carbon content in ABG understory

Enter the values obtained in this process into the calculation worksheet (in Excel) to estimate Aboveground Understory carbon content expressed in Mg C/ha.



Fig.11b Complete weight record



Fig.10.Weight of tray plus dry sample

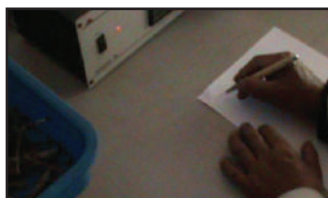


Fig.11a. Record weight of dry sample

Table 1. Estimating carbon content of plant materials( shrubs, herbs and litter) using the template in Microsoft excel

| Field work   |                                       |   |                                       | Lab work  |  |                                      | Calculated shrub biomass |                             | I. Calculated shrub C |
|--|---------------------------------------|---|---------------------------------------|---|--|--------------------------------------|--------------------------|-----------------------------|-----------------------|
| A. All shrubs + tared total wet wt<br>(kg / 25m <sup>2</sup> ) | B. Tared wt (eg, ropes, tarp)<br>(kg) | C. All shrub wet wt<br>(kg / 25m <sup>2</sup> ) | D. Subsample of shrubs<br>wet wt (kg) | E. ~ 70C oven dried subsample wt<br>dry wt (kg) | F. dry wt / wet wt fraction<br>(dry kg / wet kg) | G. Field Moisture Content<br>H2O (%) | H. Estimated dry wt      |                             |                       |
| (Enter actual measurement here)                                | (Enter actual measurement here)       | formula: A-B                                    | (Enter actual measurement)            | (Enter actual measurement)                      | formula: E/D                                     | formula: ((D-E)/D)*100               | formula: F*C             | formula: H(i)/25*10000/1000 | formula: H(ii)/2      |

Table 2. Estimating carbon content of herbs sampled from 1meter by 1meter plot

| Field work   |                                       |   |                                      | Lab work  |  |                                      | Calculated shrub biomass |                            | I. Calculated herb C |
|--|---------------------------------------|---|--------------------------------------|---|--|--------------------------------------|--------------------------|----------------------------|----------------------|
| A. All herbs + tared total wet wt<br>(kg / 1m <sup>2</sup> ) | B. Tared wt (eg, ropes, tarp)<br>(kg) | C. All herb wet wt<br>(kg / 1m <sup>2</sup> ) | D. Subsample of herbs<br>wet wt (kg) | E. ~ 70C oven dried subsample wt<br>dry wt (kg) | F. dry wt / wet wt fraction<br>(dry kg / wet kg) | G. Field Moisture Content<br>H2O (%) | H. Estimated dry wt      |                            |                      |
| (Enter actual measurement here)                              | (Enter actual measurement here)       | formula: A-B                                  | (Enter actual measurement)           | (Enter actual measurement)                      | formula: E/D                                     | formula: ((D-E)/D)*100               | formula: F*C             | formula: H(i)/1*10000/1000 | formula: H(ii)/2     |

Table 3. Estimating carbon content of organic matter/detritus of forest floor(FF) from 20cm by 20cm plot

| Field work  |  |  |  | Lab work  |  |  | Calculated shrub biomass   |  | I. Calculated herb C |
|---|--|--|--|---|--|--|----------------------------|--|----------------------|
| A. All FF + tared total wet wt<br><br>(g/400cm <sup>2</sup> ) | B. Tared wt (eg, ropes, tarp)<br><br>(g) | C. All shrub wet wt<br><br>(g/400cm <sup>2</sup> ) | D. Subsample of shrubs<br><br>wet wt (g) | E. ~ 70C ovendried subsample wt<br><br>dry wt (g) | F. dry wt / wet wt fraction<br><br>(dry g / wet g) | G. Field Moisture Content<br><br>H2O (%) | H. Estimated dry wt        |  |                      |
|   |  |  |  |   |  |  | i. (g/400cm <sup>2</sup> ) | ii. (Mg/ha)                                  | i. (Mg C/ha)         |
| (Enter actual measurement here)                               | (Enter actual measurement here)          | formula:<br>A-B                                    | (Enter actual measurement)               | (Enter actual measurement)                        | formula:<br>E/D                                    | formula:<br>((D-E)/D)*100                | formula:<br>F*C            | formula:<br>H(i)/400*10000*<br>10000/1000000 | formula:<br>H(ii)/2  |



# Chapter III:

## Procedures for soil organic carbon analysis

### 1. Soil sample sorting by the client

1.1 Upon reaching soil samples in the Laboratory Fig 12, a client should first do the sorting of the samples Fig 13 and arrange them in sequence with their identification code or labels as listed in the Sample Submission Form Fig 14. If not already filled in the Sample Submission Form the client has to fill-in the FORM as he or she wants the results to appear. If needed SPAL technician can provide help and guide how samples to be arranged.

1.2 Once the samples are arranged the client can then hand over them to the Laboratory Sample Receptionist.



Fig.12. Bulk soil samples



Fig.13. Sorting samples

### 2. Samples reception and registration by Laboratory staff

2.1 The Laboratory Sample Receptionist should check whether they are properly packed in a plastic bag and clearly labeled with identification code. Any sample not properly labeled, or labeled but not clear will be rejected unless the person handing in could clarify the confusion.

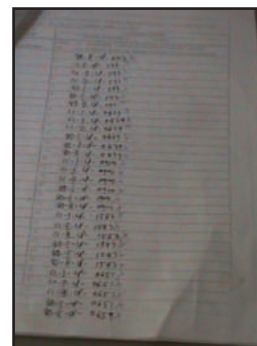


Fig.14 Sample Submission Form



Fig.15. Soil Registration Log book

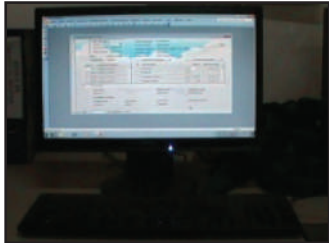


Fig.16. SPAL's Database computer



Fig.17. Samples for drying



Fig.18. Drying soil samples

- 2.2 The Sample Receptionist will also cross check whether the number of samples and their identification Code listed in the Sample Submission Form tally with the physical count and actual code on the samples. Any mismatch should be instantly corrected.
- 2.3 The Sample receptionist should also check if the contact address provided in the Form is correct.
- 2.4 The sample receiving staff should then transcript the information from the Sample Submission Form Fig 14 into the Sample Registration Form Fig 15. This in turn should be entered into the SPAL's Database computer Fig 16 to assign Lab Number and Job Identification Number (JID). This JID may be used by the client for future reference for any inquiry. Once laboratory number is assigned to respective samples the samples are then processed as in section 3 below.
- 2.5 The completed Soil Submission Form is signed and photocopied. One copy is filed and the other is returned to the originator. This is to facilitate the better communication

### 3. Soil samples processing

- 3.1 Transfer the samples Fig 17 into the drying room and dry the samples as in Fig 18.
- 3.2 Once air drying is completed Segregate the rocks greater than 2 mm from the soil sample Fig 19.
- 3.3 Weigh the stones + tarewt (eg, a tray or a plastic bag)Record the weight of the stones as in Fig 19.
- 3.4 Take the separate weight of the container.





## 4. Soil moisture content determination and correction factor

### 4.1 Principle

The moisture content is determined by drying the sample overnight at 105 °C. The moisture correction factor is calculated from the moisture content.

### 4.2. Appartus

Analytical balance with accuracy 0.001 g  
Drying Oven  
Desiccators filled with Silica gel  
Labeled Aluminium tins with lids

### 4.3. Procedure

- 4.3.1. Dry an Aluminium tin with its corresponding lid at 105 °C for 1 hour
- 4.3.2. Cool in a desiccators for 45 minutes
- 4.3.3. weigh the tin + lid with 0.001 g accuracy (Weight W g)
- 4.3.4. Transfer approximately 5.0 g of wet soil sample into the tin Fig 26.
- 4.3.5. Weigh the tin + lid + wet sample with 0.001 g accuracy (Weight A g)
- 4.3.6. Remove the lid from the tin.
- 4.3.7. Dry the tin with the sample overnight at 105 °C Fig 27.

4.3.8. Next morning remove the tin from the oven and close it with the corresponding lid Fig 28.

4.3.9. Cool down in a desiccators for about 45 minutes or until a room temperature is attained (see 5.1) Fig 29.

4.3.10. Weigh the tin + Lid + drysample with 0.001 g accuracy (Weight B g)  
Fig 30

4.3.11. Record the readings in the Moisture Logbook Fig 31.

4.3.12. Enter the recorded values obtained into the SPAL Database Fig 32 for calculating soil carbon results on dry basis.

### 4.4. Calculations

Percentage moisture air-dry basis

$$\frac{(A-B) \times 100}{(A-W)}$$

Where:

A = Weight of the moisture tin + lid + wet sample

B = Weight of the moisture tin + lid + dried sample

W = Weight of the moisture tin + lid

Record the results in 3 decimal places.

Moisture correction factor (Mc)

$$\frac{100}{(100 - Mc)}$$

Where:

M = moisture percentage as calculated above.  
Record the results in 3 decimals

#### 4.5. Remarks

As oven dried soil attracts atmospheric moisture very rapidly, the tins once removed from the oven should be placed immediately in the desiccators, and should not be left in the desiccators for more than 45 minutes.

The moisture correction factor is used to calculate results on oven-dry basis when soil samples are not dried at 105 °C before analysis. The results of the analyses are multiplied with the correction factor to obtain results on oven-dry basis.

#### 4.6. Reference: USDA, SCS, Analysis No.3A1



Fig.26 Weighing for moisture

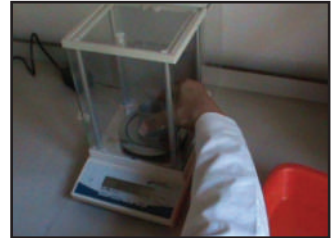


Fig.30. Weighing dry sample



Fig.27. Oven drying for moisture



Fig.31. Dry weight recorded



Fig.28. Oven dried sample

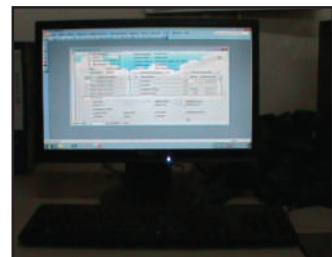


Fig.32. Moisture calculation in Database.



Fig.29. Cooling in desiccators

## 5. Analysis for soil Organic Carbon using Walkley and Black Method

Soil organic carbon can either be analyzed by using WALKLEY AND BLACK Method or with CHN Analyzer. First the procedure for Walkley and Black method is described in detail and then a brief description for CHN628 Analyzer method is presented for soil carbon analysis in section 6.

### 5.1. Principle

Soil organic matter is oxidized at a temperature of approximately 120 °C with a mixture of potassium dichromate and concentrated sulphuric acid (wet combustion). The excess potassium dichromate is titrated against ammonium ferrous sulphate with diphenylamine as an indicator. Before titration, phosphoric acid is added to form a complex with the interfering iron(III), providing a sharper colour change of the indicator.

### 5.2. Apparatus

Burette, graduated in 0.01 ml  
Erlenmeyer flask of 500 ml  
Mortar and pestle  
Sieve of 2.00 mm  
Dispenser  
Measuring cylinder  
Volumetric flask of 1000.0 ml  
Wooden pad or a white tile

## 5.3. Reagents

5.3.1. Potassium dichromate 0.1667 M (1.000 N)

1. Dry about 55 g of  $K_2Cr_2O_7$  (AR grade) at 150-200 °C for 2 hours.
2. Cool in a desiccators.
3. Weigh 49.033 g of dried  $K_2Cr_2O_7$
4. Transfer into a 1000.0 ml volumetric flask
5. Dissolve in distilled water.
6. Make up to the mark and mix.

5.3.2. Concentrated Sulphuric acid  $H_2SO_4$  with a concentration not less than 95%

5.3.3 Concentrated O-Phosphoric acid  $H_3PO_4$  with a concentration of 85%

5.3.4 Ammonium Ferrous Sulphate 0.5 M

- 1 Weigh 196 g of  $(NH_4)_2SO_4 \cdot FeSO_4 \cdot 6H_2O$
- 2 Transfer it into a 1000.0 ml volumetric flask.
- 3 Add about 700 ml of distilled water and dissolve
- 4 carefully add 20 ml of concentrated sulphuric acid and swirl
- 5 Let it cool
- 6 Make up to the mark with water and mix  
Store the solution in a brown bottle (See 5.6.1)  
For each series of

determinations standardise against  $K_2Cr_2O_7$  (See procedure below)

#### 5.3.5 Diphenylamine indicator 1%

1. Pour 20 ml of distilled water into a 250 ml Erlenmeyer flask
- 2 Carefully add 50 ml of concentrated sulphuric acid to the Erlenmeyer flask
- 3 Swirl and let it cool
- 4 Carefully add another 50 ml of concentrated sulphuric acid to the Erlenmeyer flask
- 5 Swirl and let it cool
- 6 Add 1.0 g of Diphenylamine
- 7 Dissolve and mix Store the solution in a brown bottle (See 5.6.2)

#### 5.4 Procedure

A series consists of 14 samples, including one blank and one standard sample

5.4.1 Grind about 2 g of sample in a mortar so that everything passes through a 0.5 mm sieve.

5.4.2 Weigh 1.000 g of the sieved sample into a 500 ml Erlenmeyer flask (accuracy 0.001 g) Fig 33.

5.4.3 Add 10.0 ml of 0.1667 M potassium dichromate solution with a burette or pipette and swirl the Erlenmeyer flask gently to disperse the soil in the solution.

5.4.4 Start with the blank. Add carefully 20 ml of concentrated sulphuric

acid with a dispenser or measuring cylinder Fig 34.

5.4.5 Swirl for 1 minute

5.4.6 Leave the Erlenmeyer flask on a wooden pad for 30 minutes

5.4.7 Add 200 ml of distilled water with a measuring cylinder.

5.4.8 Add 10 ml of phosphoric acid with a dispenser or a measuring cylinder.

5.4.9 Add 1.0 ml of diphenylamine indicator

5.4.10 Titrate the sample with ammonium ferrous sulphate. Near the end point the colour changes to deep violet-blue. Slow down the titration by adding the ammonium ferrous sulphate drop wise Fig 35. At the end point the colour changes sharply to brilliant green (See 5.6.3), and record the reading Fig.36 and use for calculating the carbon analysis in SPAL Database Fig 37.

5.4.11 Only for the Blanks:

After the titration of the blank, record the burette reading.

Add again 10.0 ml of 0.1667 M  $K_2Cr_2O_7$  to the titrated blank.

Titrate again with ammonium ferrous sulphate.

This second titration is for standardising the ammonium ferrous sulphate solution.





Fig.33. Weighing for carbon analysis



Fig.34. Adding sulphuric acid to sample

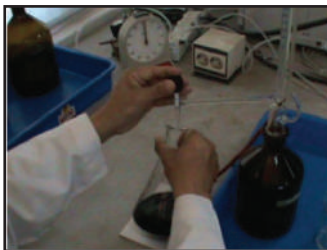


Fig.35. Titration

|    |       |        |  |
|----|-------|--------|--|
| 1  | 5488  | 15.00  |  |
| 2  | 54876 | 12.90  |  |
| 3  | 54897 | 4.60   |  |
| 4  | 54898 | 8.90   |  |
| 5  | 54899 | 13.40  |  |
| 6  | 54900 | 2.20   |  |
| 7  | 54901 | 11.80  |  |
| 8  | 54902 | Repeat |  |
| 9  | 54903 | CR     |  |
| 10 | 54904 | 8.10   |  |
| 11 | 54905 | 12.20  |  |
| 12 | 54906 | 15.70  |  |
| 13 | 54907 | 18.40  |  |

Fig.36. Record the reading

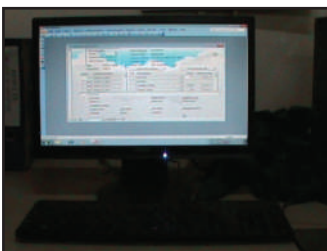


Fig.37. Result calculation in Database

## 5.5 Calculations

Normality of the ammonium ferrous solution:

$$N = \frac{F \times K}{T}$$

Where:

F = ml of potassium dichromate added to the titrated blank

K = Normality of the potassium dichromate

T = ml of ammonium ferrous sulphate used for the second titration of the blank

N = Normality of the ammonium ferrous sulphate.

Record results in 4 decimal places.

Organic carbon content:

$$\% O.C = \frac{(B-A) \times N \times 0.396 \times Mc}{W}$$

Where:

B = ml ammonium ferrous sulphate used for the first blank titration.

A = ml ammonium ferrous sulphate used for the sample.

N = Normality of the ammonium ferrous sulphate.

W = Weight of the sample in g

Mc = Moisture correction factor

0.396 is a constant factor (See 5.6.5)

Record results in 1 decimal place.

5.5.1. Enter the soil carbon values obtained in this process to estimate soil carbon content expressed in Mg C/ha using the Excel calculation worksheet shown in Table 4.

**Table 4. Estimation of C/ha for mineral soil from 10cm by 5cm from depth of 10 cm. ( The C/ha will be estimated for three soil depth: 0-10cm, 10-20cm and 20-30cm)**

| Field work                         |   |                                 |                                 |  |                  |   |   | Lab work                                      |                                    |                                  |                  |                                 |  |                                   |
|------------------------------------|---|---------------------------------|---------------------------------|--|------------------|---|---|---|------------------------------------|----------------------------------|------------------|---------------------------------|--|-----------------------------------|
| A                                  | B   | C                               | D                               | E                                      | F                | G   | H   | I   | J                                  | K                                | L                | M                               | N.(i) N(ii)                                  |                                   |
| >2mm rock etc + tared total wet wt | >2mm rock sample tared wt (eg, bucket, etc) | >2mm rock wt                    | <2mm soil + tared total wet wt  | <2mm sample tared wt (eg, bucket, etc) | <2mm soil wet wt | <2mm soil subsample for Carbon analysis                   | <2mm soil subsample for moisture calc.                      | <2mm soil subsample ~ 105C oven dried wt      | <2mm soil dry wt / wet wt fraction | <2mm soil Field Moisture Content | <2mm soil dry wt | <2mm soil Lab Carbon analysis   | Calculated Soil (<2mm) Organic C             |                                   |
| (kg/10x5x 10cm)                    | (kg)  | (kg/10x5x 10cm)                 | (kg/10x5x 10cm)                 | (kg)                                   | (kg/10x5x 10cm)  | wet wt(g)   | wet wt(g)   | dry wt(g)                                     | (dry g/ wet g)                     | H2O (%)                          | (kg/10x5x 10cm)  | Soil Organic C (%)              | i. (kg C / 50cm <sup>2</sup> x 0-10cm depth) | ii. (Mg C/ ha for 0-10cm depth)   |
| Mini                               | (Enter actual measurement here)             | (Enter actual measurement here) | (Enter actual measurement here) | (Enter actual measurement here)        | formula: D-E     | (Enter actual measurement of sub sample taken for carbon) | (Enter actual measurement of sub sample taken for moisture) | (Enter actual measurement of dry weight here) | formula: I/H                       | formula: (H-I)/H*100             | formula: F*J     | (Enter the actual SOC measured) | formula: M/100*L                             | formula: N/50* 10000* 10000/ 1000 |

## 5.6 Remarks

5.6.1 Instead of ammonium ferrous sulphate, a 0.5 M ferrous sulphate solution may be used. Ferrous sulphate however, is less stable than ammonium sulphate.

- weigh 139 g of  $\text{FeSO}_4 \cdot 7\text{H}_2\text{O}$  into a 1.0 Erlenmeyer flask
- Add about 700 ml of distilled water and dissolve.
- Carefully add 15 ml of concentrated sulphuric acid and mix
- Transfer the solution into a 1000.0 ml volumetric flask.
- make up to the mark with distilled water and mix.

Store the solution in a brown bottle.

5.6.2 An alternative indicator is barium diphenylamine sulphonate (0.16 g in 100 ml of distilled water). This indicator is easier to prepare.

5.6.3 If a sample uses less than 4 ml of 0.5 M ammonium ferrous sulphate, corresponding to a Carbon content of 2.5% or higher, repeat the determination using only 0.500 g of soil sample.

5.6.4 If a sample with 0.500 g still uses less than 4 ml of 0.5 M ammonium ferrous sulphate, corresponding to Organic Carbon content of 5.0% or higher, repeat the determination using only 0.500 g of soil sample and double the amounts of potassium dichromate, sulphuric acid, phosphoric acid and distilled water.

5.6.5 The factor 0.396 is derived from:

- Incomplete combustion (1.32)
- Equivalent weight of carbon (3)
- Conversion from milliequivalents to equivalents (1000)

d Conversion to percentage  
(100)

$$\text{i.e., } 0.396 = \frac{1.32 \times 3 \times 100}{1000}$$

#### 5.6.6 References

Black et.al. Part 2 page 1372-1375  
USDA, SCS Analysis No. 6A1

### 6. Analysis of Soil Organic Carbon using CHN628 Analyzer

The CHN628 Series Elemental Determinator is used to determine nitrogen, carbon, and hydrogen in organic matrices. The instrument utilizes a combustion technique using IR detectors for C (CO<sub>2</sub>), H (H<sub>2</sub>O) and S (SO<sub>2</sub>) and a TCD for N. The instrument features custom Microsoft Windows-based software operated through an external PC to control the system operation and data management. The standard operating procedure for the equipment is as follows:

1. Turn on the CHN628 Analyzer Fig.38
2. Turn on the computer and open 628 Series Software Fig 39
3. Open all gases (oxygen, helium, and compressed air) at gas cylinders Fig.40
4. Wait for the CHN628 furnace and reduction tube temperature to reach to its set

temperature 950 °C, and 850 °C.

5. Leave the instrument to be stabilised about 1-2 hours after the CHN628 furnace and reduction tube reach to its temperature set point.
6. Run the system check at Diagnostic page and ensure all leak check is performed.



Fig.38. CHN628 Analyzer

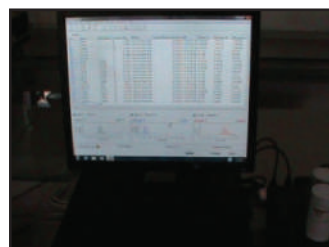


Fig.39. PC to control the system

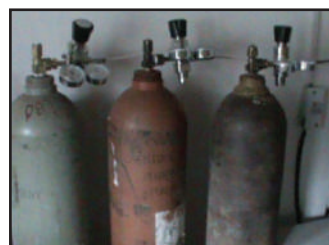


Fig.40 Gases cylinders





**Table 5. Estimation of C/ha for mineral soil from 10cm by 5cm from depth of 10 cm. ( The C/ha will be estimated for three soil depth: 0-10cm, 10-20cm and 20-30cm)**

| Field work                         |   |                                 |                                 |  |                  |   |   | Lab work                                      |                                    |                                  |                  |                                 |                                  |  |
|------------------------------------|---|---------------------------------|---------------------------------|--|------------------|---|---|---|------------------------------------|----------------------------------|------------------|---------------------------------|----------------------------------|--|
| A                                  | B   | C                               | D                               | E                                      | F                | G   | H   | I   | J                                  | K                                | L                | M                               | N.(i)                            | N(ii)                                    |
| >2mm rock etc + tared total wet wt | >2mm rock sample tared wt (eg, bucket, etc) | >2mm rock wt                    | <2mm soil + tared total wet wt  | <2mm sample tared wt (eg, bucket, etc) | <2mm soil wet wt | <2mm soil subsample for Carbon analysis                   | <2mm soil subsample for moisture calc.                      | <2mm soil subsample ~ 105C oven dried wt      | <2mm soil dry wt / wet wt fraction | <2mm soil Field Moisture Content | <2mm soil dry wt | <2mm soil Lab Carbon analysis   | Calculated Soil (<2mm) Organic C |  |
| (kg/10x5x 10cm)                    | (kg)  | (kg/10x5x 10cm)                 | (kg/10x5x 10cm)                 | (kg)                                   | (kg/10x5x 10cm)  | wet wt(g)   | wet wt(g)   | dry wt(g)                                     | (dry g/ wet g)                     | H2O (%)                          | (kg/10x5x 10cm)  | Soil Organic C (%)              | i. (kg C / 50cm2 x 0-10cm depth) | ii. (Mg C/ ha for 0-10cm depth)          |
| Mini                               | (Enter actual measurement here)             | (Enter actual measurement here) | (Enter actual measurement here) | (Enter actual measurement here)        | formula: D-E     | (Enter actual measurement of sub sample taken for carbon) | (Enter actual measurement of sub sample taken for moisture) | (Enter actual measurement of dry weight here) | formula: I/H                       | formula: (H-I)/H*100             | formula: F*J     | (Enter the actual SOC measured) | formula: M/100*L                 | formula: N/50)* 10000* 10000/ 1000/ 1000 |

## 6.17 Reference

CHN628 Carbon/Hydrogen/Nitrogen Determinator Instruction Manual; Version 1.3x, Part Number 200-731.